


EUROIMMUNUS

EUROLINE Antigen Profiles

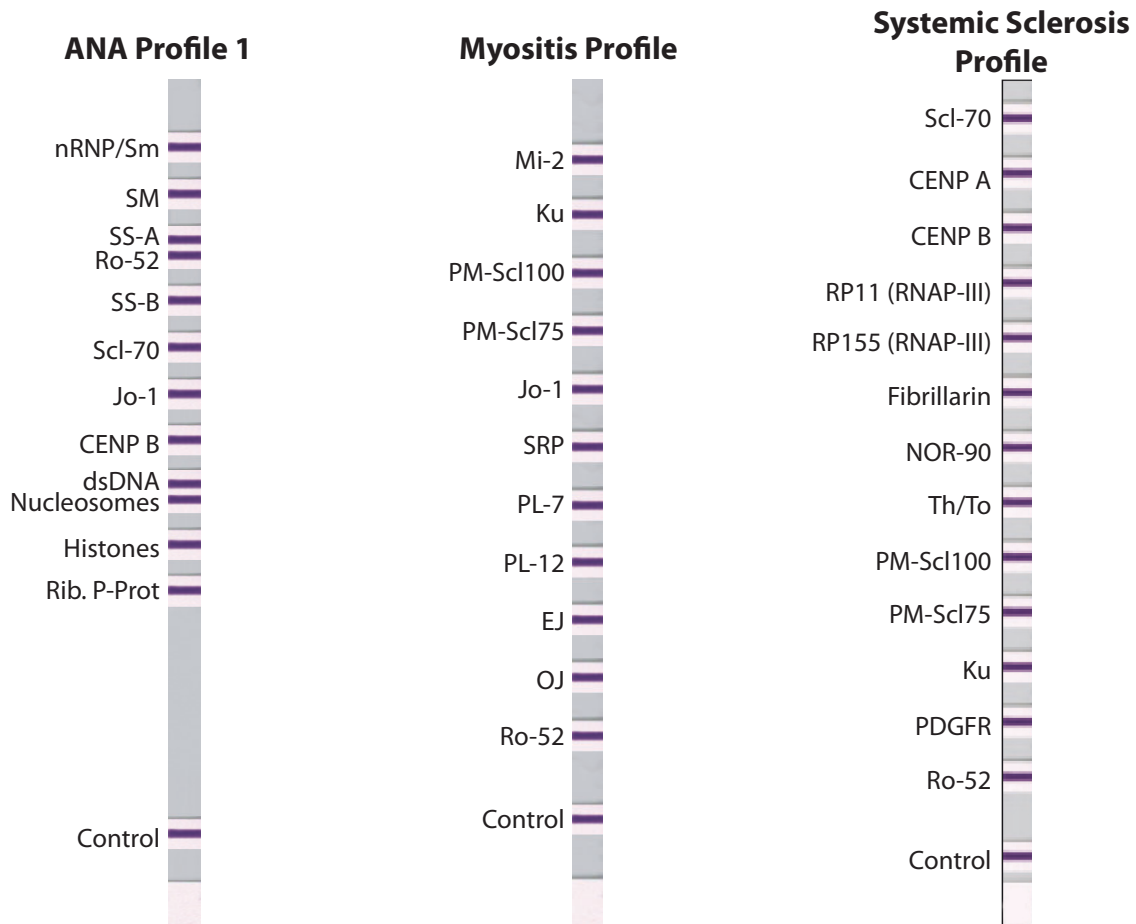

Technical Data

Antigen substrates: Membrane Chips printed with individual lines of purified, biochemically characterized antigens. Each group of antigens is coated onto a separate membrane chip, providing optimal efficiency of antibody detection for each antigen.

Sample dilution: Serum or plasma; 1:101 in sample buffer

Test procedure: 30 min/30 min/10 min (at room temperature)

Automation options: Compatible with all commercial blot processing systems.

EUROLineScan Program available for computer-based evaluation of results

Kit format: 16 membrane strips and all necessary reagents

